

THE DENTAL  
SOLUTIONS  
COMPANY™



Intego

Performance you can rely on

[dentsplysirona.com](http://dentsplysirona.com)





## Your reliable partner

Intego is the best treatment center of its class. It is based on the innovative power and award-winning design of our existing treatment centers. As a reliable partner for an efficient practice, Intego has all the features that provide optimal support: ergonomic design for better workflows, high quality standards in processing and design, and a great value – this is why Intego is our performance class.

### Great value

Excellent performance and innovative functions at a great price, customized options for great value.

### Ergonomic design

Designed for intuitive handling and optimized, integrated workflows.

### Quality made in Germany

Defined quality assurance processes, durable materials, and sturdy, future-proof design.

### Intego TS (OTP arm design)

- Ergonomic positioning of the dentist element due to a wide range support arm system
- Functional tray system with detachable tray that can be positioned flexibly
- Integrated, comprehensive infection prevention concept for optimal results



# Great value you can rely on

Your work day is complex and places many different demands on the treatment center you choose. Intego gives you the perfect solution: Choose from two treatment concepts that you can configure flexibly. Whether you select a whip arm or a hanging hoses design - every option offers optimal performance at a fair price.

### Your benefits at a glance

- Treatment center with a compact water unit for optimal use of space.
- EasyPad user interface for easy, intuitive operation.
- Compact assistant element with good freedom of movement.
- Controlled water quality with a self-contained fresh water supply bottle.
- Available as the TS model with hanging hoses or CS with whip arms.

### CS - The whip arm design



The **whip arm design** is ideal for the widest range of treatment situations with four-handed or two-handed dentistry. The dentist element can be rotated by up to 240° for optimal positioning. The maximum distance between the dentist element and instruments is 90 cm.



**The elegant world**

Classic, timeless combinations for an exceptional high-end look.




**The natural world**

Soft, warm shades for a soothing treatment atmosphere.



**The vital world**

Vivid colors and bright tones for a modern, lively atmosphere.

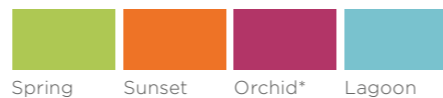

\* Available colors for Lounge upholstery.

# Ergonomic design you can rely on

The treatment center is the heart of your practice. Intego perfectly combines appearance and functionality – so you can work comfortably under optimal conditions. In addition to their harmonious design, coordinated components and variety of color combinations, Intego treatment centers meet the strict criteria of our ergonomics concept. Intego fits the style of your practice – and your taste.

**The upholstery options**



**Premium upholstery**

The high-quality, seamless upholstery ensures comfort even during long treatments.



**Lounge upholstery**

Soft, pleasant surface material and comfortable padding for utmost relaxation. Contoured surface for perfect positioning.

**Award-winning design**

Dentsply Sirona has earned prestigious design prizes for many of its products, including the red dot design award, the iF Product Design Award Gold, the iF communication design award, and the Good Design Award.



## Perfect results

Every treatment is a complex procedure. Our ergonomic products provide optimal access to the patient and smooth treatment workflows with a high level of patient comfort. Intego combines every dimension of ergonomics to make treatment easier for you. Achieve perfect results even at the end of a long working day - with your expertise and the support of our treatment center.



### 1 Intuitive sitting

- **More freedom of movement**  
Thin backrest for ideal freedom of leg movement, flexible positioning of the foot control
- **Optimized instrument placement**  
Perfect coordination of dentist and assistant elements, flexible positioning of the additional tray in the whip arm and hanging hose models
- **Comfortable working stools**  
Hugo, Carl, Paul

### 2 Comfortable positioning

- **Easy adjustment of the treatment center**  
Flexible height adjustment, swiveling cuspidor
- **Greater patient comfort**  
Gentle chair lift drive, armrests, Vario footrest, comfortable upholstery, ErgoMotion

### 3 Optimum visibility

- **Perfect illumination of the treatment area**  
LEDview Plus or LEDlight operating light, LED instrument light
- **Flexible headrests**  
Flat or double articulating headrest
- **Perfect intraoral images**  
SiroCam F, AF or AF+, 22" HD monitor on the treatment center

### 4 Integrated workflow

- **Intuitive operation**  
EasyPad user interface or EasyTouch touchscreen
- **Integrated tabletop device functions**  
Integrated torque control for endodontics, ApexLocator on the EasyTouch
- **Modern equipment**  
Innovative instruments, modern patient communication system Sivision, USB, SD card and Ethernet interface

## Perfect results

- Clinically perfect
- Esthetically perfect
- Perfect for the success of your practice

- 1 Sprayvit E**  
Optimal visibility due to integrated LED light.
- 2 BL ISO E motor**  
Sterilizable brushless motor with speed and torque control during treatment. Especially lightweight, compact and with a wide speed range.
- 3 SiroSonic L scaler**  
Handpiece and tips can be sterilized. Controlled, straight oscillation with electronic power and water control.
- 4 SiroCam F or AF or AF+**  
Detailed intraoral images with the easy to operate SiroCam F or AF/AF+ with autofocus.

## Intuitive operation

Difficult treatments require your full concentration. This is why Intego is intuitive – with two state-of-the-art, intuitive user interfaces and innovative instruments for stress-free treatment. So you can concentrate on what’s really important – your patient and a perfect treatment result.

### State-of-the-art, intuitive interfaces



**EasyPad**  
Clear and intuitive user interface – standard in Intego.




**EasyTouch**  
This user interface can be optionally integrated. More operating comfort thanks to intuitive touchscreen technology and additional operating and display options.

### Integration possibilities for endodontic treatments



**ApexLocator**  
An ApexLocator can also be integrated into the EasyTouch to allow for permanent control of the working depth during endodontic treatment and a reduction of X-ray exposure.

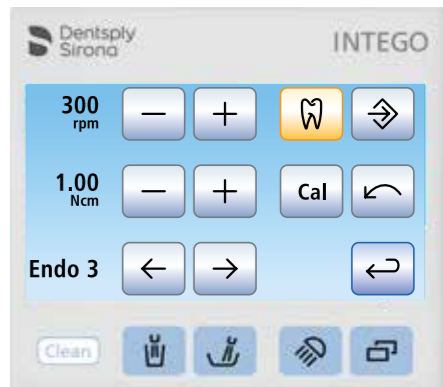

**Integrated Torque Control**  
Speed and torque are adjustable, and the motor stops automatically when the set maximum torque is reached. This ensures greater safety and perfect workflow for endodontic treatment.



### Reciprocal file Systems

Intego integrates the Dentsply Mailefer and VDW reciprocal file systems, enabling a fast and easy use of common file systems.

- WaveOne®
- WaveOne® Gold  
Reciprocating Files
- Reciproc®  
One file endo



## Impressive patient comfort

Your patients will also appreciate this comfort: The Intego patient chair has an ergonomic design and, together with the headrest, ensures that your patients are especially comfortable. The chair lowers far enough for even shorter patients to can get in and out easily. And there is no doubt – patients who are comfortable feel more secure and remain loyal to your practice.



**NEW:** Dentsply Sirona's unique ErgoMotion for the Intego family

### Your benefits

- Better positioning for maxillar or mandibular treatment
- Less overstretching of patient's head
- No overstretching and compressing of patient's spine
- No readjustment of headrest necessary

### Double articulating headrest

Quickly adjustable to various treatment situations as well as patients with varying statures.

### Premium/Lounge upholstery

The high-quality upholstery ensures comfort even during long treatments.

### **NEW:** ErgoMotion

The tilting of the backrest matches the movement of the patient's spine to prevent stretching and compression effects.

### Armrests

Your patients can rest their arms comfortably.

### Lifting system

The unique chair lift drive gently moves patients into the proper position.

### Vario footrest

The Vario footrest that comes with the Premium upholstery can be easily extended to provide even very tall patients with ergonomic positioning.

## Safety you can rely on

Dentists all over the world have the highest demands for the safety of their patients and practice team. The Intego infection prevention concept is designed so you can comply quickly and easily with all infection control regulations. The assistant element and water unit designs make infection prevention procedures even easier and less time consuming. All of this guarantees perfect infection control.



**Cuspidor**  
Removable and swiveling,  
easy to clean,  
and easily accessible.

**Purge function**  
Automatically purges  
all instrument waterlines  
on the dentist and  
assistant side.

**Suction systems**  
Compatible with  
all conventional  
wet and dry suction  
systems, including  
amalgam separator.

### Intego

The Compact water unit and assistant element can be adapted to many different demands and requirements by freely selectable suction systems and infection prevention options.



**Fresh water bottle**  
For controlled water quality. Can  
also be filled with Dentosept P  
solution for permanent disinfection  
and manual sanitization.



**Automatic suction hose cleaning**  
Cleaning adaptor for automatic  
flushing of saliva ejector and spray  
aspirator.



**Smooth surfaces**  
The largely gap-free design of the  
surfaces ensures maximum infection  
control in the shortest time possible.



## Flexible into the future

In a good practice, customized settings add to your comfort. This is why you can configure the Intego exactly the way you need for optimal support of your treatment workflow. But flexibility means even more to us – and you can prepare Intego with technical options and interfaces so you are ready for future applications.

### Many models

#### Choice of two headrests

The flat headrest can be adjusted quickly and flexibly. Optionally, you can select a double articulating headrest that is easier to position and allows you optimal visibility even in difficult-to-see areas.



Flat headrest



Double articulating headrest

#### Two operating light options – always optimal illumination

Two choices of LED lights for optimal illumination of the operating field. Our premium operating light LEDview Plus offers with its polychrome LEDs and sophisticated technical features natural light in perfection.



LEDlight



LEDview Plus

#### Assistant element: Comfort or Compact

Both assistant elements are very space-saving and ensure maximum freedom of movement. With the Comfort model, you can work comfortably even without an assistant.




Comfort assistant element



Compact assistant element

#### Ultramodern equipment

The SiroCam F/AF/AF+ autofocus cameras combine superior image quality with easy operation.

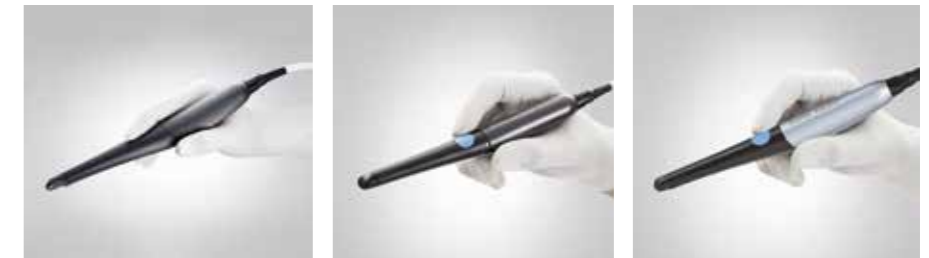

	SiroCam F	SiroCam AF	SiroCam AF+
<b>Focus</b>	Fixed Focus	Auto Focus	Auto Focus
<b>Image format 16:10</b>	✓	✓	✓
<b>Pixels (usable)</b>	0.42 megapixel	0.65 megapixel	1.00 megapixel
<b>Workspace</b>	✓ 5–20 mm	✓ 2–130 mm	✓ 2–130 mm
<b>Resolution live-image</b>	824 x 514	824 x 514	824 x 514
<b>Resolution still-image</b>	824 x 514	1024 x 640	1276 x 796

### Future-proof design

#### Interfaces: flexible into the future

A network interface in the patient chair and a USB port in the dentist or assistant element allow you to integrate new functions and network your practice – troubleshooting by your service technician is also much easier.



Integrated USB interface

#### Sivision 22" AC monitor

With this HD monitor, you can explain intra-oral camera and X-ray images, software and planning views, videos and PowerPoint presentations to your patient directly at the treatment center.



Sivision 22" AC monitor

# Quality you can rely on

We don't want to just satisfy our customers, we want to impress them. This is why quality, performance, and service are key for us. Dentsply Sirona treatment units combine high quality with real innovations that perceptibly enhance treatment comfort and workflow.

## Reliability and durability



Precise, elaborate quality tests and robust materials make Intego a reliable treatment center with an above-average service life.

## High quality



Intego immediately shows our attention to detail. The materials used - robust, shapeable plastics, solid steel and aluminum components on the inside, high-quality glass and upholstery - meet the high quality demands of a premium product.

## High standards



Dentsply Sirona is a technology leader with more than 130 years of innovation history. Our engineers and scientists in Research & Development not only pay close attention to innovation, but also quality.

## Future-proof design



At the Dentsply Sirona innovation center, future trends are anticipated and products are prepared to be equipped for the technology demands of the future.



# Individual configuration

You can configure your Intego treatment center to fit your needs and meet your own desired great value.

## Intego

Basic unit	Version	User packages	Option packages
<p><b>Patient chair</b></p> <ul style="list-style-type: none"> <li>• Premium upholstery</li> <li>• Flat headrest</li> <li>• Pneumatic instrument foot control</li> </ul> <p><b>Dentist element</b></p> <ul style="list-style-type: none"> <li>• 5 instrument positions</li> <li>• EasyPad user interface</li> <li>• Supply for 2 drill drives</li> </ul> <p><b>Assistant element</b></p> <ul style="list-style-type: none"> <li>• Compact version on a swivel arm directly mounted to the cuspidor</li> <li>• 4 handpiece positions</li> <li>• Saliva ejector and HVE</li> </ul> <p><b>Water unit</b></p> <ul style="list-style-type: none"> <li>• Compact design</li> <li>• Standard wet suction</li> <li>• Ready for the integration of a fresh water bottle</li> <li>• Cuspidor with manual swivel function</li> </ul>	<p><b>TS version with hanging hoses</b></p>  <ul style="list-style-type: none"> <li>• The height of the dentist element can be adjusted and held in place with a pneumatic brake</li> <li>• Ergonomic design with flexibly positionable dentist element and movable tray holder</li> <li>• Ample instrument area on the dentist element and tray holder</li> </ul> <p><b>CS version with whip arms</b></p>  <ul style="list-style-type: none"> <li>• Ideal for the widest range of treatment situations, regardless of two-handed or four-handed treatment</li> <li>• Optimum reach of instrument hoses: up to 90 cm</li> <li>• The height of the dentist element can be adjusted and held in place with a pneumatic brake</li> </ul>	<p><b>Start</b></p> <ul style="list-style-type: none"> <li>• Standard 3-way syringe</li> <li>• Light source in the dentist element</li> <li>• 2 turbine hoses</li> <li>• LEDlight*</li> </ul> <p><b>Innovation</b></p> <ul style="list-style-type: none"> <li>• C+ electronic foot control</li> <li>• Sprayvit E in the dentist element</li> <li>• Light source in the dentist element</li> <li>• 1 turbine hose</li> <li>• 1 BL ISO E micro-motor with supply electronics</li> <li>• Sirosonic L scaler</li> <li>• LEDview Plus</li> </ul> <p><b>Ergo</b></p> <ul style="list-style-type: none"> <li>• ErgoMotion</li> <li>• Chair lift system 185 kg</li> <li>• 2 armrests</li> <li>• Double articulated headrest</li> </ul> <p><b>Endo Pro**</b></p> <ul style="list-style-type: none"> <li>• Integrated torque control</li> <li>• ApexLocator</li> <li>• Endo 6:1 contra-angle handpiece</li> <li>• e-Control</li> </ul> <p><b>Advance</b></p> <ul style="list-style-type: none"> <li>• C+ electronic foot control</li> <li>• Sprayvit E in the dentist element</li> <li>• Light source in the dentist element</li> <li>• 1 turbine hose</li> <li>• 1 BL ISO E micro-motor with supply electronics</li> <li>• LEDlight*</li> </ul> <p><b>Touch Plus</b></p> <ul style="list-style-type: none"> <li>• C+ electronic foot control</li> <li>• 1 BL ISO E micro-motor with supply electronics</li> <li>• 1 turbine hose</li> <li>• Light source in the dentist element</li> <li>• EasyTouch user interface</li> <li>• LEDview Plus</li> </ul> <p><b>** LEDview Plus at additional charge</b>  <b>** Only with user interface EasyTouch</b></p>	<p><b>Sivision Ready</b></p> <ul style="list-style-type: none"> <li>• USB interface for the dentist or assistant element</li> <li>• Network interface</li> </ul> <p><b>Sivision Start</b></p> <p><b>with/without monitor</b></p> <ul style="list-style-type: none"> <li>• USB interface for the dentist or assistant element</li> <li>• Network interface</li> <li>• SiroCam F (SiroCam AF or AF+ at additional charge)</li> <li>• Monitor mount on the light support column</li> <li>• Optional: 22" AC monitor</li> </ul>



# Your individual options

Here you can view the complete list of individual options, and a checklist of what is included with basic units.

Patient Chair	
2 treatment positions, S and O position	■
Integrated junction box	■
Chair lift system 140 kg	■
Flat headrest	■
Double articulating headrest	<input type="checkbox"/> Er
Premium upholstery	■
Lounge upholstery	<input type="checkbox"/>
ErgoMotion	<input type="checkbox"/> Er
4-way foot switch	<input type="checkbox"/>
Booster cushion for small patients	<input type="checkbox"/>
Right armrest	<input type="checkbox"/> Er
Left armrest	<input type="checkbox"/> Er
Valve set for external devices	<input type="checkbox"/>
Vacuum pressure reduction kit	<input type="checkbox"/>
Cover installation kit	<input type="checkbox"/>
Separate junction box	<input type="checkbox"/>
Adapter plate	<input type="checkbox"/>

Dentist element	
Height-adjustable support arm with pneumatic brake	■
5 instrument positions	■
Hydro-pneumatic supply for 2 drill drives w/o fiber optics	■
Pneumatic instrument foot control	■
C+ electronic foot control	<input type="checkbox"/> A, I, T
EasyPad user interface with 2 user-specific settings	■
EasyTouch user interface with 4 user-specific settings	<input type="checkbox"/> T
Sprayvit E with light	<input type="checkbox"/> A, I
Sprayvit E with light (water button on the left)	<input type="checkbox"/> A, I
Standard 3-way syringe	<input type="checkbox"/> S
Turbine hose with light	<input type="checkbox"/> 2xS, A, I, T
Borden turbine hose	<input type="checkbox"/>
1st BL ISO E motor with hose	<input type="checkbox"/> A, I, T
2nd BL ISO E motor with hose	<input type="checkbox"/>
SiroSonic L scaler	<input type="checkbox"/> I
Third drill drive	<input type="checkbox"/>
Intra-oral light source kit for turbine and Sprayvit E	<input type="checkbox"/> S, A, I, T
Torque control for endodontic treatments	<input type="checkbox"/> E
eControl: Speed range of 100 - 40,000 rpm	<input type="checkbox"/> E
Integrated ApexLocator	<input type="checkbox"/> E
Endo 6:1 contra-angle handpiece	<input type="checkbox"/> E
Intra-oral X-ray viewer	<input type="checkbox"/>
Purge and manual sanitization function	<input type="checkbox"/>
Silicone pad for TS dentist element	<input type="checkbox"/>
Additional CS tray holder for 1 standard size tray	<input type="checkbox"/>
Additional TS/CS tray holder for 2 standard size trays	<input type="checkbox"/>
Silicone pad for CS dentist element	<input type="checkbox"/>
Removable silicone sleeves for handles	<input type="checkbox"/>

Assistant element	
Compact assistant element, swivel mounted directly on the cuspidor arm	■
Comfort assistant element, ideal for two-handed dentistry	<input type="checkbox"/>
4 instrument positions	■
Saliva ejector and HVE	■
Sprayvit E with light	<input type="checkbox"/>
Sprayvit E with light (water button on the left)	<input type="checkbox"/>
Standard 3-way syringe	<input type="checkbox"/>
2nd HVE kit	<input type="checkbox"/>
Satelec Mini LED curing light	<input type="checkbox"/>
Touch pad	<input type="checkbox"/>
Light source for Sprayvit E	<input type="checkbox"/>
Silicone pad for Comfort assistant element	<input type="checkbox"/>

Water unit	
Prepared for fresh water bottle	■
Cuspidor with manual swivel function	■
Water unit without cuspidor	<input type="checkbox"/>
Prepared for standard wet suction	■
Prepared for standard wet suction w/automatic suction control	<input type="checkbox"/>
Multi-station valve	<input type="checkbox"/>
Cuspidor valve for standard wet suction	<input type="checkbox"/>
Automatic separator	<input type="checkbox"/>
Amalgam separator	<input type="checkbox"/>
Fresh water bottle w/o switch to public water supply	<input type="checkbox"/>
Fresh water bottle with switch to public water supply	<input type="checkbox"/>
Water heater	<input type="checkbox"/>
Water flush suction hose cleaning system	<input type="checkbox"/>

Operating light options	
LEDlight	<input type="checkbox"/> S, A
LEDview Plus	<input type="checkbox"/> I, T
LEDview Plus ceiling light	<input type="checkbox"/>
Adaptation for LEDview Plus/Heliident Plus ceiling combination	<input type="checkbox"/>
Adaptation for Heliident Plus unit model	<input type="checkbox"/>

Patient communication	
USB interface in the dentist or assistant element	<input type="checkbox"/> SiR, SiS
SiroCam F	<input type="checkbox"/> SiS
SiroCam AF	<input type="checkbox"/>
SiroCam AF+	<input type="checkbox"/>
Monitor mount on the light support column	<input type="checkbox"/> SiS
22" AC monitor (with integrated power supply unit)*	<input type="checkbox"/> SiS*
Network interface	<input type="checkbox"/> SiR, SiS

- Included in basic unit
- Option
- Not possible/applicable
- SiR Sivation Ready package
- SiS Sivation Start package
- SiS\* Sivation Start package with monitor
- S Start package
- A Advance package
- I Innovation package
- E Endo Pro package
- Er Ergo package
- T Touch Plus package



Hugo working stool



Carl working stool



Paul working stool (pictured with additional foot ring)



Additional swiveling tray holder for CS dentist element



Additional swiveling tray holder for TS dentist element



Double articulating headrest



Satelec Mini LED curing light



X-ray viewer



Additional holder for camera (for TS dentist element)



Comfort assistant element



C+ electronic foot control



Flat headrest



Valve set for external devices



Adaptation for Heliident Plus unit model



22" AC monitor



Intego without cuspidor



LEDlight



LEDview Plus



SiroCam AF



SiroCam AF+

## Dentsply Sirona

Sirona Dental Systems GmbH  
Fabrikstraße 31, 64625 Bensheim, Deutschland  
[dentsplysirona.com](http://dentsplysirona.com)

## Procedural Solutions

Preventive  
Restorative  
Orthodontics  
Endodontics  
Implants  
Prosthetics

## Enabling Technologies

CAD/CAM  
Imaging  
Treatment Centers  
Instruments